

ENG

Instructions for Use:
**HUMAN sINTERLEUKIN-2
RECEPTOR ELISA**

Catalogue number:
RAF106R

For research use only!

 **BioVendor**
R&D[®]



BioVendor – Laboratorní medicína a.s.

Karásek 1767/1, 621 00 Brno, Czech Republic

+420 549 124 185

info@biovendor.com

sales@biovendor.com

www.biovendor.com

1. INTENDED USE	3
2. SUMMARY	3
3. PRINCIPLES OF THE TEST	4
4. REAGENTS PROVIDED	5
5. STORAGE INSTRUCTIONS – ELISA KIT	5
6. SPECIMEN COLLECTION AND STORAGE INSTRUCTIONS	5
7. MATERIAL REQUIRED BUT NOT SUPPLIED	5
8. PRECAUTIONS FOR USE	6
9. PREPARATION OF REAGENTS	7
10. TEST PROTOCOL	9
11. CALCULATION OF RESULTS	11
12. LIMITATIONS	12
13. PERFORMANCE CHARACTERISTICS	13
14. REAGENT PREPARATION SUMMARY	16
15. EXPLANATION OF SYMBOLS	17
16. ASSAY PROCEDURE - SUMMARY	18

HISTORY OF CHANGES

Previous version	Current version
ENG.006.A	ENG.007.A
History of changes added.	
A symbol indicating the manufacturer added.	

1. INTENDED USE

The Human sIL-2R ELISA is an enzyme-linked immunosorbent assay for the quantitative detection of Human sIL-2R. **The Human sIL-2R ELISA is for research use only. Not for diagnostic or therapeutic procedures.**

2. SUMMARY

Interleukin-2 (IL-2) has been described as a factor that promotes the growth and proliferation of Human T-cells thus occupying a pivotal role in the generation of the immune response. This proliferation of T lymphocytes is triggered by the interaction of IL-2 with its specific cell surface receptor following T lymphocyte activation. The receptor for IL-2 is composed of at least three distinct polypeptide subunits, the IL-2R, IL-2R, and the IL-2R chains giving rise to a 55-65 kD membrane bound protein.


The genes encoding IL-2 and these three receptor subunits have been cloned and their complete primary structures have been deduced. Evidence has accumulated that suggests a critical role of IL-2R in IL-2 signal transduction. IL-2R coupled tyrosine kinases are believed to play a crucial role in that signal transduction.

3. PRINCIPLES OF THE TEST

An anti-Human sIL-2R coating antibody is adsorbed onto microwells.

Figure 1


Coated Microwell

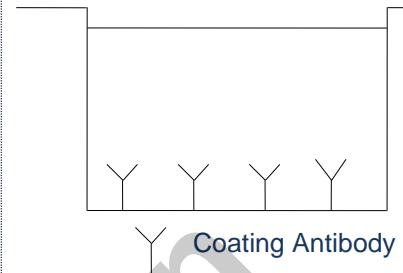

Human sIL-2R present in the sample or standard binds to antibodies adsorbed to the microwells. A HRP-conjugated anti-Human sIL-2R antibody is added and binds to Human sIL-2R captured by the first antibody.

Figure 2


First Incubation

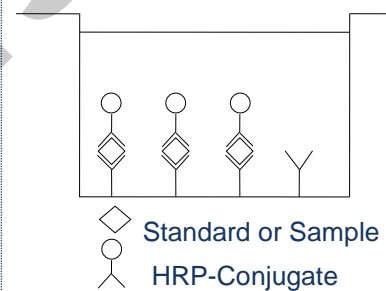

Following incubation unbound HRP-conjugated anti-Human sIL-2R is removed during a wash step, and substrate solution reactive with HRP is added to the wells.


Figure 3

Second Incubation

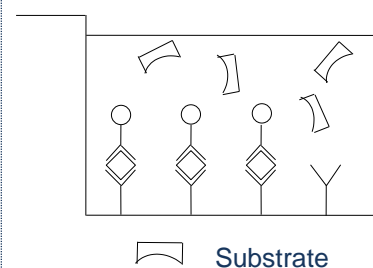

A coloured product is formed in proportion to the amount of Human sIL-2R present in the sample or standard. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 Human sIL-2R standard dilutions and Human sIL-2R concentration determined.

Figure 4

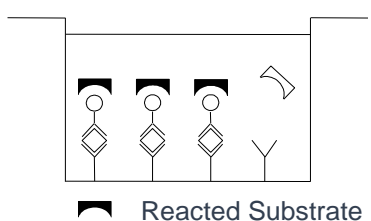

4. REAGENTS PROVIDED

- 1 aluminium pouch with a **Antibody Coated Microtiter Strips** with monoclonal antibody to Human sIL-2R
- 1 vial (70 µl) **HRP-Conjugate** anti-Human sIL-2R monoclonal antibody
- 2 vials Human sIL-2R **Standard** lyophilized, 40 ng/ml upon reconstitution
- 1 vial (12 ml) **Sample Diluent**
- 1 vial (5 ml) **Assay Buffer Concentrate** 20x (PBS with 1% Tween 20 and 10% BSA)
- 1 bottle (50 ml) **Wash Buffer Concentrate** 20x (PBS with 1% Tween 20)
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 2 **Adhesive Films**

5. STORAGE INSTRUCTIONS – ELISA KIT

Store kit reagents between 2°C and 8°C. Immediately after use remaining reagents should be returned to cold storage (2°C to 8°C). Expiry of the kit and reagents is stated on labels.

Expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

6. SPECIMEN COLLECTION AND STORAGE INSTRUCTIONS

Cell culture supernatant, serum and plasma (EDTA, citrate, heparin) were tested with this assay. Other biological samples might be suitable for use in the assay. Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Pay attention to a possible “**Hook Effect**” due to high sample concentrations (see chapter Calculation of results **Chyba! Nenalezen zdroj odkazů.**).

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive Human sIL-2R. If samples are to be run within 24 hours, they may be stored at 2°C to 8°C (for sample stability refer to 13.5).

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

7. MATERIAL REQUIRED BUT NOT SUPPLIED

- 5 ml and 10 ml graduated pipettes
- 5 µl to 1000 µl adjustable single channel micropipettes with disposable tips
- 50 µl to 300 µl adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform regression analysis

8. PRECAUTIONS FOR USE

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.
- Exposure to acid inactivates the conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

9. PREPARATION OF REAGENTS

Buffer Concentrates should be brought to room temperature and should be diluted before starting the test procedure.

If crystals have formed in the **Buffer Concentrates**, warm them gently until they have completely dissolved.

9.1 Wash Buffer (1x)

1. Pour entire contents (50 ml) of the **Wash Buffer Concentrate (20x)** into a clean 1000 ml graduated cylinder. Bring to final volume of 1000 ml with glass-distilled or deionized water. Mix gently to avoid foaming.
2. Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer (1x) is stable for 30 days.
3. Wash Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Wash Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

9.2 Assay Buffer (1x)

1. Pour the entire contents (5 ml) of the **Assay Buffer Concentrate (20x)** into a clean 100 ml graduated cylinder. Bring to final volume of 100 ml with distilled water. Mix gently to avoid foaming.
2. Store at 2°C to 8°C. Please note that the Assay Buffer (1x) is stable for 30 days.
3. Assay Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Assay Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

9.3 HRP-Conjugate

Please note that the HRP-Conjugate should be used within 30 minutes after dilution.

Make a 1:100 dilution of the concentrated HRP-Conjugate solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	HRP-Conjugate (ml)	Assay Buffer (1x) (ml)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

9.4 Human sIL-2R Standard


1. Reconstitute Human sIL-2R standard by addition of distilled water.
2. Reconstitution volume is stated in the Quality Control Sheet. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 40 ng/ml).
3. Allow the standard to reconstitute for 10-30 minutes. Mix well prior to making dilutions.
4. After usage remaining standard cannot be stored and has to be discarded.
5. Standard dilutions can be prepared directly on the microwell plate or alternatively in tubes (see 9.4.1).

9.4.1 External Standard Dilution

1. Label 7 tubes, one for each standard point. S1, S2, S3, S4, S5, S6, S7
2. Prepare 1:2 serial dilutions for the standard curve as follows:
Pipette 225 μ l of Sample Diluent into each tube.
3. Pipette 225 μ l of reconstituted standard (concentration = 40 ng/ml) into the first tube, labelled S1, and mix (concentration of standard 1 = 20 ng/ml).
4. Pipette 225 μ l of this dilution into the second tube, labelled S2, and mix thoroughly before the next transfer.
5. Repeat serial dilutions 5 more times thus creating the points of the standard curve (see Figure 5).

Sample Diluent serves as blank.

Figure 5 Dilute standards - tubes

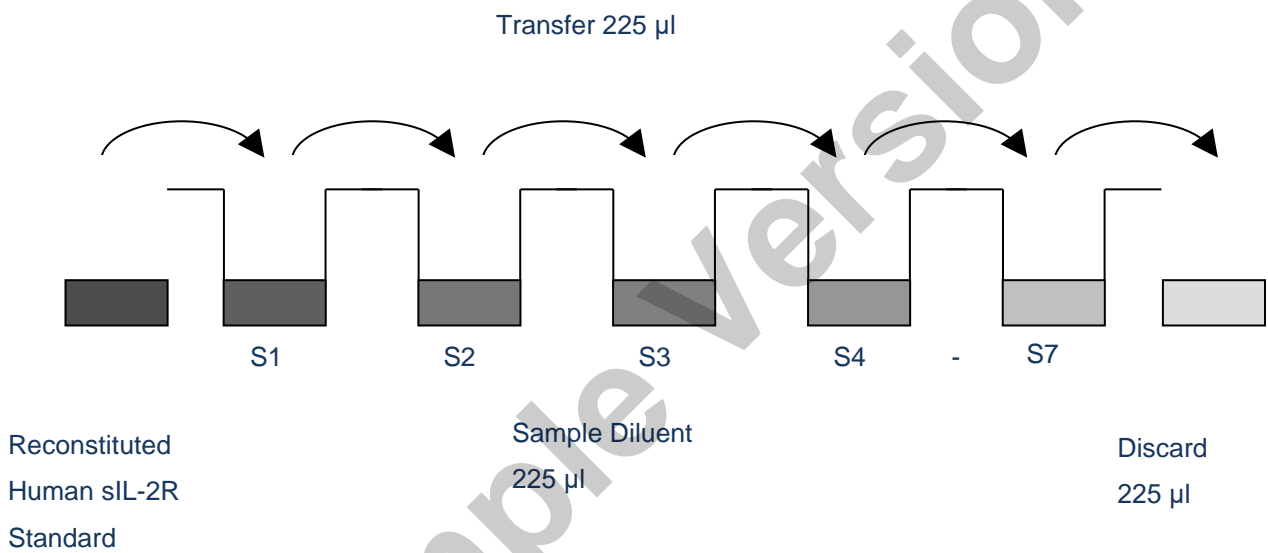

10. TEST PROTOCOL

- Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control sample should be assayed in duplicate. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2°-8°C sealed tightly.
- Wash the microwell strips twice with approximately 400 µl **Wash Buffer** per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about **10 – 15 seconds** before aspiration. Take care not to scratch the surface of the microwells. After the last wash step, empty wells and tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing. Alternatively microwell strips can be placed upside down on a wet absorbent paper for not longer than 15 minutes. **Do not allow wells to dry.**
- Standard dilution on the microwell plate** (Alternatively the standard dilution can be prepared in tubes – see 9.4.1): Add 100 µl of Sample Diluent in duplicate to all **standard wells**. Pipette 100 µl of prepared **standard** (see Preparation of Standard 9.4, concentration = 40.00 ng/ml) in duplicate into well A1 and A2 (see Table 1). Mix the contents of wells A1 and A2 by repeated aspiration and ejection (concentration of standard 1, S1 = 20.00 ng/ml), and transfer 100 µl to wells B1 and B2, respectively (see Figure 6). Take care not to scratch the inner surface of the microwells. Continue this procedure 5 times, creating two rows of Human sIL-2R standard dilutions ranging from 20.00 to 0.31 ng/ml. Discard 100 µl of the contents from the last microwells (G1, G2) used.

Figure 6 Dilute standards – microwell plate

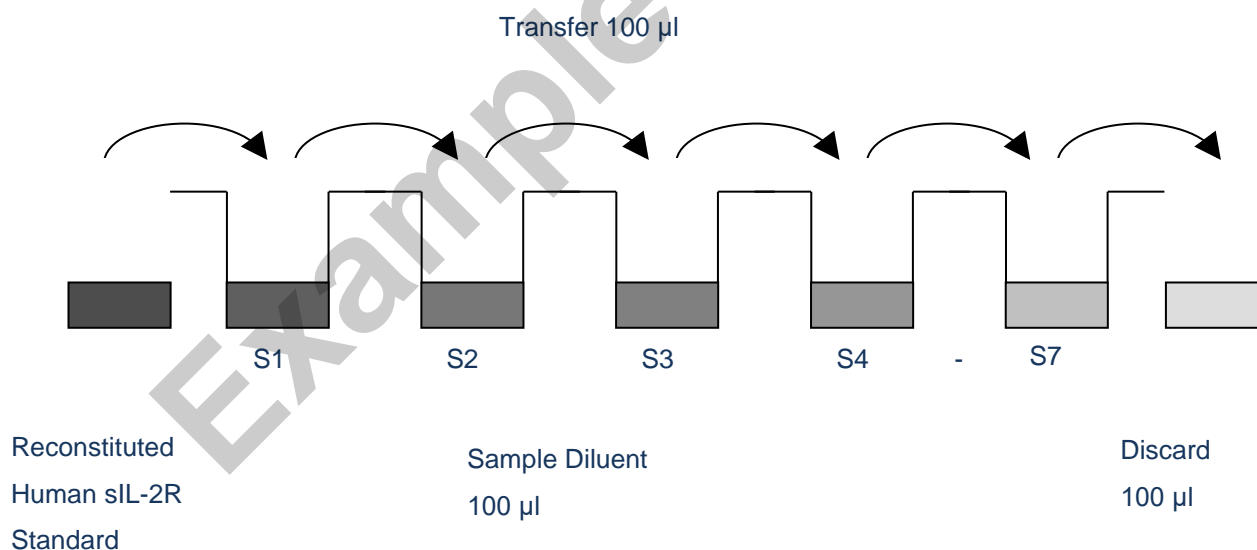


Table 1 Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	1	2	3	4
A	Standard 1 (20.00 ng/ml)	Standard 1 (20.00 ng/ml)	Sample 1	Sample 1
B	Standard 2 (10.00 ng/ml)	Standard 2 (10.00 ng/ml)	Sample 2	Sample 2
C	Standard 3 (5.00 ng/ml)	Standard 3 (5.00 ng/ml)	Sample 3	Sample 3
D	Standard 4 (2.50 ng/ml)	Standard 4 (2.50 ng/ml)	Sample 4	Sample 4
E	Standard 5 (1.25 ng/ml)	Standard 5 (1.25 ng/ml)	Sample 5	Sample 5
F	Standard 6 (0.63 ng/ml)	Standard 6 (0.63 ng/ml)	Sample 6	Sample 6
G	Standard 7 (0.31 ng/ml)	Standard 7 (0.31 ng/ml)	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

In case of an external standard dilution (see “External standard dilution“), pipette 100 µL of these standard dilutions (S1 - S7) in the standard wells according to Table 1


4. Add 100 µl of **Sample Diluent** in duplicate to the **blank wells**.
5. Add 50 µl of **Sample Diluent** to the **sample wells**.
6. Add 50 µl of each **sample** in duplicate to the **sample wells**.
7. Prepare **HRP-Conjugate** (see Preparation of HRP-Conjugate).
8. Add 50 µl of **HRP-Conjugate** to all wells.
9. Cover with an adhesive film and incubate at room temperature (18 to 25°C) for 3 hours, if available on a microplate shaker set at 400 rpm.
10. Remove adhesive film and empty wells. **Wash** microwell strips 3 times according to point 2. of the test protocol. Proceed immediately to the next step.
11. Pipette 100 µl of **TMB Substrate Solution** to all wells.
12. Incubate the microwell strips at room temperature (18° to 25°C) for about 10 min. Avoid direct exposure to intense light. The colour development on the plate should be monitored and the substrate reaction stopped (see next point of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay. It is recommended to add the stop solution when the highest standard has developed a dark blue colour. Alternatively the colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.9 – 0.95.
13. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
14. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the standards.

Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

11. CALCULATION OF RESULTS

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean value.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the Human sIL-2R concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating Human sIL-2R for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding Human sIL-2R concentration.
- If instructions in this protocol have been followed samples have been diluted 1:2 (50 μ l sample + 50 μ l Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor ($\times 2$).
- Calculation of samples with a concentration exceeding standard 1 will result in incorrect, low Human sIL-2R levels (Hook Effect). Such samples require further external predilution according to expected Human sIL-2R values with Sample Diluent in order to precisely quantitate the actual Human sIL-2R level.
- It is suggested that each testing facility establishes a control sample of known Human sIL-2R concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 7. This curve cannot be used to derive test results. Each laboratory must prepare a standard curve for each group of microwell strips assayed.

Figure 7 Representative standard curve for Human sIL-2R ELISA. Human sIL-2R was diluted in serial 2-fold steps in Sample Diluent.

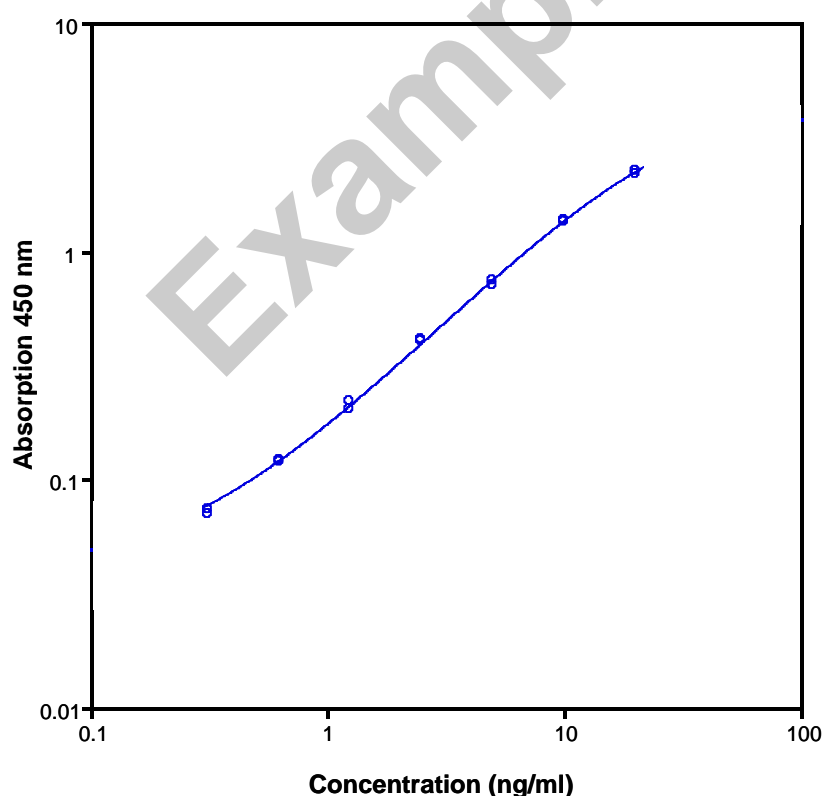

Note: Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

Table 1 Typical data using the Human sIL-2R ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Human sIL-2R Concentration (ng/ml)	O.D. at 450 nm	Mean O.D. at 450 nm	C.V.(%)
1	20.00	2.239	2.216	3.25
		2.193		
2	10.00	1.366	1.360	0.85
		1.354		
3	5.00	0.710	0.729	2.62
		0.747		
4	2.50	0.403	0.407	0.49
		0.410		
5	1.25	0.204	0.212	1.06
		0.219		
6	0.63	0.120	0.121	0.14
		0.122		
7	0.31	0.071	0.073	0.21
		0.074		
Blank	0	0.010	0.011	0.14
		0.012		

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.

12. LIMITATIONS

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of radioimmunotherapy has significantly increased the number of patients with Human anti-mouse IgG antibodies (HAMA). HAMA may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the sample.

13. PERFORMANCE CHARACTERISTICS

13.1 Sensitivity

The limit of detection of Human sIL-2R defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 0.27 ng/ml (mean of 6 independent assays).

13.2 Reproducibility

13.2.1 Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 6 serum samples containing different concentrations of Human sIL-2R. 2 standard curves were run on each plate. Data below show the mean Human sIL-2R concentration and the coefficient of variation for each sample (see Table 2). The calculated overall intra-assay coefficient of variation was 7.2%.

Table 2 The mean Human sIL-2R concentration and the coefficient of variation for each sample

Sample	Experiment	Mean Human sIL-2R Concentration (ng/ml)	Coefficient of Variation (%)
1	1	2.60	11
	2	2.53	7
	3	2.01	3
2	1	1.47	10
	2	1.54	7
	3	1.32	5
3	1	1.47	9
	2	1.56	7
	3	1.30	9
4	1	1.72	2
	2	1.80	7
	3	1.64	7
5	1	1.43	8
	2	1.52	11
	3	1.33	4
6	1	3.05	10
	2	3.42	5
	3	2.67	4

13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 6 serum samples containing different concentrations of Human sIL-2R. 2 standard curves were run on each plate. Data below show the mean Human sIL-2R concentration and the coefficient of variation calculated on 18 determinations of each sample (see Table 3). The calculated overall inter-assay coefficient of variation was 9.8%.

Table 3 The mean Human sIL-2R concentration and the coefficient of variation of each sample

Sample	Mean Human sIL-2R Concentration (ng/ml)	Coefficient of Variation (%)
1	2.381	13.6
2	1.440	8.1
3	1.443	9.1
4	1.803	9.2
5	1.430	6.5
6	3.048	12.3

13.3 Spiking Recovery

The spiking recovery was evaluated by spiking 4 levels of Human sIL-2R into different pooled normal Human serum samples. Recoveries were determined in 3 independent experiments with 4 replicates each. The amount of endogenous Human sIL-2R in unspiked serum was subtracted from the spike values.

The overall mean recovery was 75%.

13.4 Dilution Linearity

4 serum samples with different levels of Human sIL-2R were analysed at serial 2 fold dilutions with 4 replicates each.

The recovery ranged from 77% to 107% with an overall recovery of 86% (see Table 4).

Table 4

Sample	Dilution	Expected Human sIL-2R Concentration (ng/ml)	Observed Human sIL-2R Concentration (ng/ml)	Recovery of Expected Human sIL-2R Concentration (%)
1	1:2	--	1.27	--
	1:4	6.37	5.92	93
	1:8	3.18	2.58	81
	1:16	1.59	1.28	81
2	1:2	--	17.0	--
	1:4	8.50	9.08	107
	1:8	4.25	3.83	90
	1:16	2.13	1.86	87
3	1:2	--	18.6	--
	1:4	9.30	8.35	90

	1:8	4.65	3.69	79
	1:16	2.33	1.78	77
4	1:2	--	16.72	--
	1:4	8.36	6.79	81
	1:8	4.18	3.49	83
	1:16	2.09	1.84	88

13.5 Sample Stability

13.5.1 Freeze-Thaw Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20°C and thawed 5 times, and the Human sIL-2R levels determined. There was no significant loss of Human sIL-2R immunoreactivity detected by freezing and thawing.

13.5.2 Storage Stability

Aliquots of serum samples (spiked or unspiked) were stored at -20°C, 2-8°C, room temperature (RT) and at 37°C, and the Human sIL-2R level determined after 24, 48 and 96 h. There was no significant loss of Human sIL-2R immunoreactivity detected during storage under above conditions.

13.6 Specificity

The interference of circulating factors of the immune system was evaluated by spiking these proteins at physiologically relevant concentrations into a Human IL-2R positive serum. There was no crossreactivity detected.

13.7 Expected Values

Panels of 40 serum as well as EDTA, citrate and heparin plasma samples from randomly selected apparently healthy donors (males and females) were tested for Human IL-2R. The levels measured may vary with the sample collection used.

Table 5

Sample Matrix	Number of Samples Evaluated	Range (ng/ml)	Mean (ng/ml)	Standard Deviation (ng/ml)
Serum	40	1.9 – 13.1	4.7	2.6
Plasma (EDTA)	40	0.9 – 8.1	2.4	1.3
Plasma (Citrate)	40	1.8 - 7.8	3.8	1.65
Plasma (Heparin)	40	2.12 – 8.0	4.0	1.7

14. REAGENT PREPARATION SUMMARY

14.1 Wash Buffer (1x)

Add **Wash Buffer Concentrate** 20x (50 ml) to 950 ml distilled water.

Number of Strips	Wash Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

14.2 Assay Buffer (1x)

Add **Assay Buffer Concentrate** 20x (5 ml) to 95 ml distilled water.

Number of Strips	Assay Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

14.3 HRP-Conjugate










Make a 1:100 dilution of **HRP-Conjugate** in Assay Buffer (1x):

Number of Strips	HRP-Conjugate (ml)	Assay Buffer (1x) (ml)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

14.4 Human sIL-2R Standard

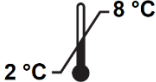Reconstitute lyophilized **Human sIL-2R standard** with distilled water. (Reconstitution volume is stated in the Quality Control Sheet.)

15. EXPLANATION OF SYMBOLS

	Catalogue number
	Batch code
	Caution
	Use by date
	Temperature limit
	Manufacturer
 <p data-bbox="258 1182 466 1205">www.biovendor.com</p>	Read electronic instructions for use - eIFU
	The content is sufficient for 96 tests
	Biological risks

16. ASSAY PROCEDURE - SUMMARY

1. Determine the number of microwell strips required.
2. Wash microwell strips twice with Wash Buffer.
3. Standard dilution on the microwell plate: Add 100 μ l Sample Diluent, in duplicate, to all standard wells. Pipette 100 μ l prepared standard into the first wells and create standard dilutions by transferring 100 μ l from well to well. Discard 100 μ l from the last wells. Alternatively external standard dilution in tubes: Pipette 100 μ l of these standard dilutions in the microwell strips.
4. Add 100 μ l Sample Diluent, in duplicate, to the blank wells.
5. Add 50 μ l Sample Diluent to sample wells.
6. Add 50 μ l sample in duplicate, to designated sample wells.
7. Prepare HRP-Conjugate.
8. Add 50 μ l HRP-Conjugate to all wells.
9. Cover microwell strips and incubate 3 hours at room temperature (18°C to 25°C).
10. Empty and wash microwell strips 3 times with Wash Buffer.
11. Add 100 μ l of TMB Substrate Solution to all wells.
12. Incubate the microwell strips for about 10 minutes at room temperature (18°C to 25°C).
13. Add 100 μ l Stop Solution to all wells.
14. Blank microwell reader and measure colour intensity at 450 nm.

Note: If instructions in this protocol have been followed samples have been diluted 1:2 (50 μ l sample + 50 μ l Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 2).



BioVendor – Laboratorní medicína a.s.

Karásek 1767/1, 621 00 Brno, Czech Republic

+420 549 124 185

info@biovendor.com

sales@biovendor.com

www.biovendor.com



FOLLOW US