

BioVendor

Research
and Diagnostic Products



Mouse Progranulin ELISA

Product Data Sheet

Cat. No.: RAG014R

For Research Use Only

CONTENTS

| | | |
|-----|-------------------------------------|----|
| 1. | INTENDED USE | 3 |
| 2. | HANDLING, STORAGE | 3 |
| 3. | INTRODUCTION | 3 |
| 4. | TEST PRINCIPLE | 4 |
| 5. | TECHNICAL HINTS | 4 |
| 6. | REAGENT SUPPLIED | 5 |
| 7. | MATERIALS REQUIRED BUT NOT SUPPLIED | 5 |
| 8. | PREPARATION OF REAGENTS | 5 |
| 9. | PREPARATION OF SAMPLES | 7 |
| 10. | ASSAY PROCEDURE | 8 |
| 11. | CALCULATIONS | 9 |
| 12. | PERFORMANCE CHARACTERISTICS | 10 |
| 13. | TROUBLESHOOTING | 12 |
| 14. | REFERENCES | 13 |

**»» This kit is manufactured by:
BioVendor – Laboratorní medicína a.s.**

»» Use only the current version of Product Data Sheet enclosed with the kit!

1. INTENDED USE

The Mouse Progranulin ELISA kit is to be used for the *in vitro* quantitative determination of mouse progranulin in serum and cell culture supernatant. This ELISA Kit is for research use only.

2. HANDLING, STORAGE

- Reagent must be stored at 2-8°C when not in use.
- Plate and reagents should be at room temperature before use.
- Do not expose reagents to temperatures greater than 25°C.

3. INTRODUCTION

The granulins are a family of cysteine-rich polypeptides, some of which have growth modulatory activity. The widespread occurrence of granulin mRNA in cells from the hematopoietic system and in epithelia implies important functions in these tissues. All 4 known human granulin-like peptides are encoded in a single precursor, progranulin, which has a highly conserved 12-cysteine backbone defining a consensus sequence that is repeated 7 times (1).

Progranulin is a 593-amino acid glycoprotein, the mRNA of which is expressed in many epithelial cells both *in vitro* and *in vivo*. He and Bateman demonstrated that overexpression of the progranulin gene in SW-13 adrenal carcinoma cells and MDCK nontransformed renal epithelia resulted in transfection-specific secretion of progranulin, acquired clonogenicity in semisolid agar, and increased mitosis in monolayer culture, whereas diminution of progranulin gene expression impaired growth of these cells (2). They proposed that the rate of growth of some epithelia is proportional to the level of intrinsic progranulin gene expression, and that elevated progranulin gene expression confers a transformed phenotype on epithelial cells including anchorage independence *in vitro* and growth as tumors in nude mice. They also found that in murine transcutaneous puncture wounds, progranulin mRNA is expressed in the inflammatory infiltrate and is highly induced in dermal fibroblasts and endothelia following injury. When applied to a cutaneous wound, progranulin increased the accumulation of neutrophils, macrophages, blood vessels, and fibroblasts in the wound. It acted directly on isolated dermal fibroblasts and endothelial cells to promote division, migration, and the formation of capillary-like tubule structures, concluding that progranulin is, therefore, probably a wound-related growth factor (3). Due to its tumor-promoting activity and multiple functions on hematopoietic cells and endothelial cells, progranulin can be a novel marker for oncology and hematology.

4. TEST PRINCIPLE

This assay is a sandwich Enzyme Linked-Immunosorbent Assay (ELISA) for quantitative determination of mouse progranulin in biological fluids. A polyclonal antibody specific for progranulin has been precoated onto the 96-well microtiter plate. Standards and samples are pipetted into the wells for binding to the coated antibody. After extensive washing to remove unbound compounds, progranulin is recognized by the addition of a biotinylated polyclonal antibody specific for progranulin (Detection Antibody). After removal of excess biotinylated antibody, HRP labeled streptavidin (STREP-HRP) is added. Following a final washing, peroxidase activity is quantified using the substrate 3,3',5,5'-tetramethylbenzidine (TMB). The intensity of the color reaction is measured at 450 nm after acidification and is directly proportional to the concentration of progranulin in the samples.

5. TECHNICAL HINTS

- It is recommended that all standards, controls and samples be run in duplicate.
- Do not combine leftover reagents with those reserved for additional wells.
- Reagents from the kit with a volume less than 100 µl should be centrifuged.
- Residual wash liquid should be drained from the wells after last wash by tapping the plate on absorbent paper.
- Crystals could appear in the 10X solution due to high salt concentration in the stock solutions. Crystals are readily dissolved at room temperature or at 37°C before dilution of the buffer solutions.
- Once reagents have been added to the 16-well strips, DO NOT let the strips DRY at any time during the assay.
- Keep Substrate Solution protected from light.
- The Stop Solution consists of phosphoric acid. Although diluted, the Stop Solution should be handled with gloves, eye protection and protective clothing.

6. REAGENT SUPPLIED

| Kit Components | Quantity |
|---|------------------|
| 1 plate coated with mouse progranulin Antibody | 6x16-well strips |
| 2 bottles Wash Buffer 10X | 30 ml |
| 2 bottle ELISA Buffer 10X | 30 ml |
| 1 vial Detection Antibody | 20 μ l |
| 1 vial HRP Labeled Streptavidin (lyophilized) | 2 μ g |
| 1 vial mouse progranulin Standard (lyophilized) | 16 ng |
| 1 bottle TMB Substrate Solution | 12 ml |
| 1 bottle Stop Solution | 12 ml |
| 2 plate sealers (plastic film) | |
| 2 silica Gel Minibags | |

7. MATERIALS REQUIRED BUT NOT SUPPLIED

- Microtiterplate reader at 450 nm
- Calibrated precision single and multi-channel pipettes. Disposable pipette tips
- Deionized water
- Microtubes or equivalent for preparing dilutions
- Disposable plastic containers for preparing working buffers
- Plate washer: automated or manual
- Glass or plastic tubes for diluting and aliquoting standard

8. PREPARATION OF REAGENTS

NOTE: Prepare just the appropriate amount of the buffers necessary for the assay.

- **Wash Buffer 10X** has to be diluted with deionized water 1:10 before use (e.g. 50 ml Wash Buffer 10X + 450 ml water) to obtain Wash Buffer 1X.
- **ELISA Buffer 10X** has to be diluted with deionized water 1:10 before use (e.g. 20 ml ELISA Buffer 10X + 180 ml water) to obtain ELISA Buffer 1X.

- **Detection Antibody (DET)** has to be diluted to 1:1000 in ELISA Buffer 1X (10 µl DET + 10 ml ELISA Buffer 1X).

NOTE: The diluted Detection Antibody is not stable and cannot be stored!

- **HRP Labeled Streptavidin (STREP-HRP)** has to be reconstituted with 100 µl of ELISA Buffer 1X.
 - After reconstitution of STREP-HRP, prepare aliquots and store them at -20°C. Avoid freeze/thaw cycles.
 - Dilute the reconstituted STREP-HRP to the working concentration by adding 8 µl in 12 ml of ELISA Buffer 1X (1:1500).

NOTE: The diluted STREP-HRP is not stable and cannot be stored!

- **Mouse progranulin Standard (STD)** has to be reconstituted with 1 ml of deionized water.
 - This reconstitution produces a stock solution of 16 ng/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes. Mix well prior to making dilutions.

NOTE: The reconstituted standard is aliquoted and stored at -20°C.

- Dilute the standard protein concentrate (STD) (16 ng/ml) in ELISA Buffer 1X. A seven-point standard curve using 2-fold serial dilutions in ELISA Buffer 1X is recommended.
- Suggested standard points are:
8 , 4 , 2 , 1 , 0.5 , 0.25 , 0.125 and 0 ng/ml.

Dilute further for the standard curve:

| To obtain | Add | Into |
|--------------------|------------------------------------|---------------------------|
| 8 ng/ml | 300 µl of progranulin (16 ng/ml) | 300 µl of ELISA Buffer 1X |
| 4 ng/ml | 300 µl of progranulin (8 ng/ml) | 300 µl of ELISA Buffer 1X |
| 2 ng/ml | 300 µl of progranulin (4 ng/ml) | 300 µl of ELISA Buffer 1X |
| 1 ng/ml | 300 µl of progranulin (2 ng/ml) | 300 µl of ELISA Buffer 1X |
| 0.5 ng/ml | 300 µl of progranulin (1 ng/ml) | 300 µl of ELISA Buffer 1X |
| 0.25 ng/ml | 300 µl of progranulin (0.5 ng/ml) | 300 µl of ELISA Buffer 1X |
| 0.125 ng/ml | 300 µl of progranulin (0.25 ng/ml) | 300 µl of ELISA Buffer 1X |
| 0 ng/ml | 300 µl of ELISA Buffer 1X | Empty tube |

9. PREPARATION OF SAMPLES

Serum : Use a serum separator tube. Let samples clot at room temperature for 30 minutes before centrifugation for 20 minutes at 1,000xg. Assay freshly prepared serum or store serum in aliquot at $\leq -20^{\circ}\text{C}$ for later use. Avoid repeated freeze/thaw cycles.

Serum or Cell Culture Supernatant have to be diluted in ELISA Buffer 1X. Samples containing visible precipitates must be clarified before use.

NOTE: *As a starting point, 1/1'000 dilution of serum is recommended! If sample values fall outside the detection range of the assay, a lower or higher dilution may be required!*


10. ASSAY PROCEDURE

| | |
|--------------------------|--|
| <input type="checkbox"/> | <p>1. Determine the number of 16-well strips needed for the assay and insert them in the frame for current use. The extra strips should be resealed in the foil pouch bag and stored at 4°C.</p> <p>NOTE: Remaining 16-well strips coated with progranulin antibody when opened can be stored at 4°C for up to 1 month.</p> |
| <input type="checkbox"/> | <p>2. Add 100 µl of the different standards into the appropriate wells in duplicate! At the same time, add 100 µl of diluted serum, plasma or cell culture supernatant samples in duplicate to the wells (see Preparation and Storage of Reagents and Preparation of Samples).</p> |
| <input type="checkbox"/> | <p>3. Cover the plate with plate sealer and incubate for 1 hour at 37°C.</p> |
| <input type="checkbox"/> | <p>4. Aspirate the coated wells and add 300 µl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of three washes. After the last wash, complete removal of liquid is essential for good performance.</p> |
| <input type="checkbox"/> | <p>5. Add 100 µl to each well of the Detection Antibody (DET) (see Preparation and Storage of Reagents).</p> |
| <input type="checkbox"/> | <p>6. Cover the plate with plate sealer and incubate for 1 hour at 37°C.</p> |
| <input type="checkbox"/> | <p>7. Aspirate the coated wells and add 300 µl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of three washes. After the last wash, complete removal of liquid is essential for good performance.</p> |
| <input type="checkbox"/> | <p>8. Add 100 µl to each well of the diluted HRP Labeled Streptavidin (STREP-HRP) (see Preparation and Storage of Reagents).</p> |
| <input type="checkbox"/> | <p>9. Cover the plate with plate sealer and incubate for 1 hour at 37°C.</p> |
| <input type="checkbox"/> | <p>10. Aspirate the coated wells and add 300 µl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of five washes. After the last wash, complete removal of liquid is essential for good performance.</p> |
| <input type="checkbox"/> | <p>11. Add 100 µl to each well of TMB Substrate Solution (TMB).</p> |
| <input type="checkbox"/> | <p>12. Allow the color reaction to develop at room temperature (RT°C) in the dark for 10 minutes.</p> |
| <input type="checkbox"/> | <p>13. Stop the reaction by adding 100 µl of Stop Solution (STOP). Tap the plate gently to ensure thorough mixing. The substrate reaction yields a blue solution that turns yellow when Stop Solution is added.</p> |
| | <p>! CAUTION: CORROSIVE SOLUTION!</p> |
| <input type="checkbox"/> | <p>14. Measure the OD at 450 nm in an ELISA reader within 30 minutes.</p> |

11. CALCULATIONS

- Average the duplicate readings for each standard, control and sample and subtract the average blank value (obtained with the 0 ng/ml point).
- Generate the standard curve by plotting the average absorbance obtained for each standard concentration on the horizontal (X) axis vs. the corresponding progranulin concentration (ng/ml) on the vertical (Y) axis (see 10. TYPICAL DATA).
- Calculate the progranulin concentrations of samples by interpolation of the regression curve formula as shown above in a form of a quadratic equation.
- If the test samples were diluted, multiply the interpolated values by the dilution factor to calculate the concentration of mouse progranulin in the samples.

The following data are obtained using the different concentrations of standard as described in this protocol:

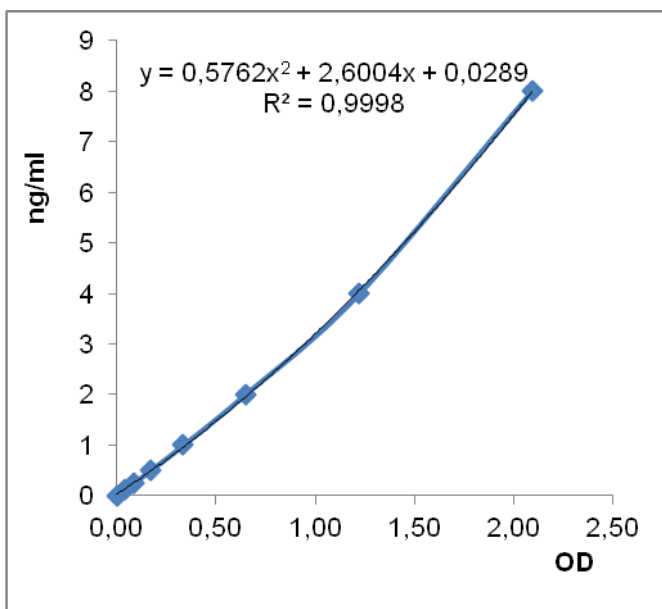

| Standard mProgranulin (ng/ml) | Optical Density (mean) |
|-------------------------------|------------------------|
| 8 | 2.09 |
| 4 | 1.22 |
| 2 | 0.65 |
| 1 | 0.33 |
| 0.5 | 0.17 |
| 0.25 | 0.09 |
| 0.125 | 0.04 |
| 0 | 0.00 |

Figure: Standard curve

12. PERFORMANCE CHARACTERISTICS

Typical analytical data of BioVendor Mouse Progranulin ELISA are presented in this chapter

- **Sensitivity** (Limit of detection)

The lowest level of progranulin that can be detected by this assay is 60 pg/ml. **NOTE:** *The Limit of detection was measured by adding two standard deviations to the mean value of 50 zero standard.*

- **Assay range**

0.125 ng/ml – 8 ng/ml

- **Specificity**

This ELISA is specific for the measurement of natural and recombinant mouse progranulin. It does not cross-react with mouse RBP4, mouse adiponectin, mouse Nampt, mouse leptin, mouse RELM- β , mouse IL-33, mouse IDO, mouse resistin, mouse clusterin, mouse ANGPTL3, mouse ANGPTL4, human vaspin, human adiponectin, human progranulin, rat Nampt.

Precision:

Intra-assay

Six samples of known concentrations of mouse progranulin were assayed in replicates 8 times to test precision within an assay.

| Samples | Means ($\mu\text{g/ml}$) | SD | CV (%) | n |
|---------|----------------------------|-------|--------|---|
| 1 | 1.089 | 0.045 | 4.10 | 8 |
| 2 | 0.938 | 0.055 | 5.83 | 8 |
| 3 | 2.289 | 0.096 | 4.18 | 8 |
| 4 | 2.096 | 0.097 | 4.62 | 8 |
| 5 | 1.689 | 0.121 | 7.19 | 8 |
| 6 | 2.459 | 0.139 | 5.67 | 8 |

Inter-assay

Six samples of known concentrations of mouse progranulin were assayed in 8 separate assays to test precision between assays.

| Samples | Means ($\mu\text{g/ml}$) | SD | CV (%) | n |
|---------|----------------------------|-------|--------|---|
| 1 | 1.154 | 0.089 | 7.75 | 8 |
| 2 | 2.351 | 0.168 | 7.13 | 8 |
| 3 | 2.037 | 0.117 | 5.77 | 8 |
| 4 | 1.693 | 0.136 | 8.03 | 8 |
| 5 | 2.210 | 0.140 | 6.35 | 8 |
| 6 | 0.951 | 0.056 | 5.93 | 8 |

Linearity

Different mouse serum samples containing progranulin were diluted several fold (1/1,000 to 1/4,000) and the measured recoveries ranged from 90% to 103%.

| Samples | Sample Dilution | Expected ($\mu\text{g/ml}$) | Observed ($\mu\text{g/ml}$) | % of Expected |
|---------|-----------------|-------------------------------|-------------------------------|---------------|
| 1 | 1 : 1,000 | 1.161 | 1.161 | 100 |
| | 1 : 2,000 | 0.581 | 0.594 | 102.29 |
| | 1 : 4,000 | 0.290 | 0.293 | 102.70 |
| 2 | 1 : 1,000 | 2.042 | 2.042 | 100 |
| | 1 : 2,000 | 1.021 | 0.951 | 93.18 |
| | 1 : 4,000 | 0.511 | 0.462 | 90.55 |
| 3 | 1 : 1,000 | 1.271 | 1.271 | 100 |
| | 1 : 2,000 | 0.636 | 0.597 | 93.95 |
| | 1 : 4,000 | 0.318 | 0.288 | 90.48 |

Recovery

When samples (serum) are spiked with known concentrations of mouse progranulin, the recovery averages 97% (range from 90% to 105%).

| Samples | Average recovery (%) | Range (%) |
|---------|----------------------|-----------|
| 1 | 97.85 | 95-105 |
| 2 | 100.38 | 95-105 |
| 3 | 92.66 | 90-100 |
| 4 | 97.68 | 95-105 |

Expected values:

Progranulin levels range in mouse samples from **1 to > 5 $\mu\text{g/ml}$** .

13. TROUBLESHOOTING

| PROBLEM | POSSIBLE CAUSES | SOLUTIONS |
|--------------------------|-------------------------------------|--|
| No signal or weak signal | Omission of key reagent | Check that all reagents have been added in the correct order. |
| | Washes too stringent | Use an automated plate washer if possible. |
| | Incubation times inadequate | Incubation times should be followed as indicated in the manual. |
| | Plate reader settings not optimal | Verify the wavelength and filter setting in the plate reader. |
| | Incorrect assay temperature | Use recommended incubation temperature. Bring substrates to room temperature before use. |
| High background | Concentration of STREP-HRP too high | Use recommended dilution factor. |
| | Inadequate washing | Ensure all wells are filling wash buffer and are aspirated completely. |
| Poor standard curve | Wells not completely aspirated | Completely aspirate wells between steps. |
| | Reagents poorly mixed | Be sure that reagents are thoroughly mixed. |
| Unexpected results | Omission of reagents | Be sure that reagents were prepared correctly and added in the correct order. |
| | Dilution error | Check pipetting technique and double-check calculations. |

14. REFERENCES


References to Progranulin:


1. Isolation and sequence of the granulin precursor cDNA from human bone marrow reveals tandem cysteine-rich granulin domains: V. Bhandari, et al.; Proc. Nat. Acad. Sci. **89**, 1715 (1992)
2. Progranulin gene expression regulates epithelial cell growth and promotes tumor growth in vivo: Z. He, et al.; Cancer Res. **59**, 3222 (1999)
3. Progranulin is a mediator of the wound response: Z. He, et al.; Nature Med. **9**, 225 (2003)

References to this product:

1. T. Philips, et al.; J. Neuropathol. Exp. Neurol. **69**, 1191 (2010)
2. T. Nedachi, et al.; Neuroscience **185**, 106 (2011)
3. O.A. Brady, et al.; Hum. Mol. Genet. **22**, 685 (2013)
4. A.J. Filiano, et al.; J. Neurosci. **33**, 5352 (2013)

NOTES

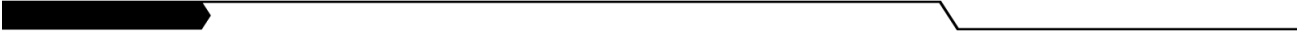

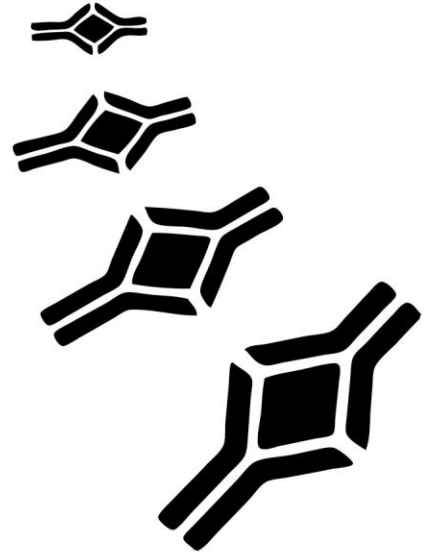

BioVendor – Laboratorní medicína a.s.

Karasek 1767/1, 621 00 Brno, Czech Republic

Phone: +420-549-124-185, Fax: +420-549-211-460

E-mail: info@biovendor.com, sales@biovendor.com

Web: www.biovendor.com

There are BioVendor branches and distributors near you.

To find the office closest to you, visit www.biovendor.com/contact

