

## PRODUCT DATASHEET

### Fibroblast Growth Factor 21 Mouse E. coli Tag free

**Cat. No.:** RD272108100-B

**Type:** Recombinant protein

**Size:** 0.1 mg

**Source:** E. coli

**Species:** Mouse

#### Description

Total 183 AA. MW: 20.1 kDa (calculated). UniProtKB acc.no. Q9JJN1. 182AA of the mouse FGF21 and one extra AA, N-terminal methionin (highlighted).

#### Other names

FGF-21, UNQ3115/PRO10196

#### Introduction to the molecule

The FGFs are a family of more than 20 small (~17–26 kDa) secreted peptides. The initial characterization of these proteins focused on their ability to stimulate fibroblast proliferation. This mitogenic activity was mediated through FGF receptors (FGFRs) 1, 2, or 3. A fourth closely related tyrosine kinase receptor (FGFR4) was able to bind the FGFs but did not lead to a mitogenic response.

FGFs modulate cellular activity via at least 5 distinct subfamilies of high-affinity FGF receptors (FGFRs): FGFR-1, -2, -3, and -4, all with intrinsic tyrosine kinase activity and, except for FGFR-4, multiple splice isoforms, and FGFR-5, which lacks an intracellular kinase domain. There is growing evidence that FGFRs can be important for regulation of glucose and lipid homeostasis. The overexpression of a dominant negative form of FGFR-1 in  $\beta$  cells leads to diabetes in mice, which thus implies that proper FGF signaling is required for normal  $\beta$  cell function and glycemia maintenance. FGFR-2 appears to be a key molecule during pancreatic development. Moreover, FGFR-4 has been implicated in cholesterol metabolism and bile acid synthesis.

FGF-19, has been shown to cause resistance to diet-induced obesity and insulin desensitization and to improve insulin, glucose, and lipid profiles in diabetic rodents. Since these effects, at least in part, are mediated through the observed changes in metabolic rates, FGF-19 can be considered as a regulator of energy expenditure.

FGF-21 is preferentially expressed in liver, but an exact knowledge of FGF-21 bioactivity and its mode of action have been lacking to date. FGF-21 is a potent activator of glucose uptake on adipocytes, protects animals from diet-induced obesity when overexpressed in transgenic mice, and lowers blood glucose and triglyceride levels when therapeutically administered to diabetic rodents.

#### Research topic

Animal studies, Diabetology - Other Relevant Products, Energy metabolism and body weight regulation

#### Amino Acid sequence

MAYPIPDSSP LLQFGGQVRQ RYLYTDDDQD TEAHLEIRED GTVVGAAHRS PESLLELKAL KPGVIQILGV KASRFLCQQP  
DGALYGSPHF DPEACSFREL LLEDGYNVYQ SEAHGLPLRL PQKDSPNQDA TSWGPRVFLP MPGLLHEPQD QAGFLPPEPP  
DVGSSDPLSM VEPLQGRSPS YAS

#### Purity

>95%

**Endotoxin**

< 0.1 EU/μg

**Formulation:**

Filtered (0,4 μm) and lyophilized in 0.5 mg/mL in 20mM TRIS, 50mM NaCl, pH 7.5

**Reconstituion:**

Add deionized water to prepare a working stock solution of approximately 0.5 mg/mL and let the lyophilized pellet dissolve completely. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture. Add DTT (0.2 mM ) and NaCl (0.1–0.15 M) before freezing to prevent potential aggregation.

**Shipping**

At ambient temperature. Upon receipt, store the product at the temperature recommended below.

**Storage, Stability/Shelf Life**

Store the lyophilized protein at –80 °C. Lyophilized protein remains stable until the expiry date when stored at –80 °C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at –80 °C for long term storage. Reconstituted protein can be stored at 4 °C for a week.

**Quality control**

BCA to determine quantity of the protein.

SDS PAGE to determine purity of the protein.


LAL to determine quantity of endotoxin.

**Applications**

ELISA, Western blotting

**Note**

This product is intended for research use only.

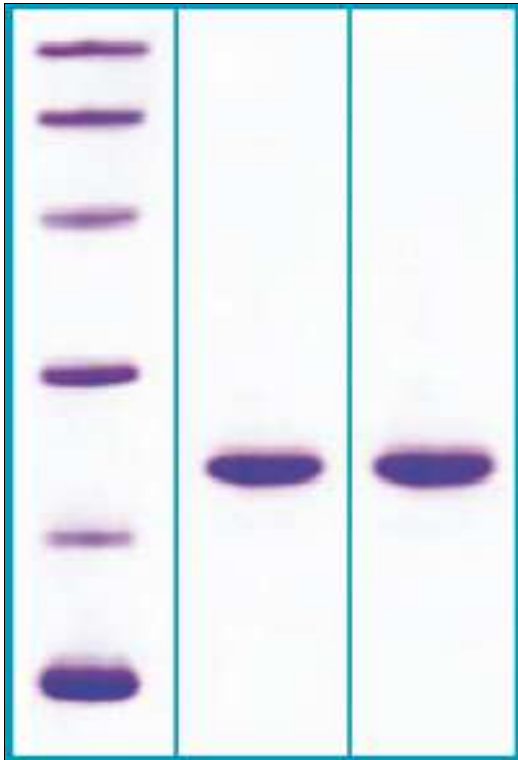

12% SDS-PAGE separation of Mouse FGF21

1. M.W. marker – 14, 21, 31, 45, 66, 97 kDa

2. reduced and heated sample, 5µg/lane

3. non-reduced and non-heated sample, 5µg/lane