

**BioVendor**

Research  
and Diagnostic Products



## HUMAN PLACENTAL LACTOGEN ELISA

Product Data Sheet

Cat. No.: RIS0013R

For Research Use Only

## CONTENTS

---

1.	INTENDED USE	3
2.	STORAGE, EXPIRATION	3
3.	INTRODUCTION	3
4.	TEST PRINCIPLE	4
5.	PRECAUTIONS	4
6.	TECHNICAL HINTS	5
7.	REAGENT SUPPLIED	6
8.	MATERIAL REQUIRED BUT NOT SUPPLIED	7
9.	PREPARATION OF REAGENTS	7
10.	PREPARATION OF SAMPLES	7
11.	ASSAY PROCEDURE	8
12.	CALCULATIONS	9
13.	PERFORMANCE CHARACTERISTICS	11
14.	LIMITATIONS OF USE	13
15.	QUALITY CONTROL	13
16.	LEGAL ASPECTS	14

**»» This kit is manufactured by:  
BioVendor – Laboratorní medicína a.s.**

**»» Use only the current version of Product Data Sheet enclosed with the kit!**

## 1. INTENDED USE

---

RIS0013R Human Placental Lactogen ELISA Immunoenzymetric assay for the research use only of human Placental Lactogen in serum. It is intended for research use only.

## 2. STORAGE, EXPIRATION

---

When stored at 2-8°C unopened reagents will retain reactivity until expiration date. Do not use reagents beyond this date.

Opened reagents must be stored at 2-8°C. Microtiter wells must be stored at 2-8°C. Once the foil bag has been opened, care should be taken to close it tightly again.

Opened kits retain activity for six weeks if stored as described above.

## 3. INTRODUCTION

---

The physiological role of human Placental Lactogen has not yet been established, but the great similarity to human growth hormone has stimulated the hypothesis about a function as a regulator of fetoplacental growth and other physiological alterations during pregnancy. It has been suggested that the maternal serum-level may reflect on "index of placental function".

Depressed levels of human Placental Lactogen are seen in association with intrauterine death, fetal distress in labor and birth asphyxia. This association is particularly strong if depressed levels are seen repeatedly, implying a chronic state of placental and therefore fetal compromise. Depressed levels are usually not present if the pregnancy has proceeded uneventfully to term.

Elevated levels of human Placental Lactogen are usually indicative of optimal pregnancy outcome in singleton pregnancies. However, high levels may indicate substantial fetal pathology in specific diseases, namely diabetes mellitus and fetal macrosomia, rhesus isoimmunization and hydrops fetalis.

## 4. TEST PRINCIPLE

---

The Biovendor human Placental Lactogen ELISA Kit is a solid phase enzyme-linked immunosorbent assay (ELISA) based on the sandwich principle.

The microtiter wells are coated with a monoclonal [mouse] antibody directed towards a unique antigenic site on an human Placental Lactogen molecule. An aliquot of patient sample containing endogenous human Placental Lactogen is incubated in the coated well with enzyme conjugate, which is a monoclonal anti- human Placental Lactogen antibody conjugated with horseradish peroxidase. After incubation the unbound conjugate is washed off.

The amount of bound peroxidase is proportional to the concentration of human Placental Lactogen in the sample.

Having added the substrate solution, the intensity of colour developed is proportional to the concentration of human Placental Lactogen in the patient sample.

## 5. PRECAUTIONS

---

### **Safety**

#### **For research use only**

All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazard in use and for disposal.

Before starting the assay, read the instructions completely and carefully. Use the valid version of the package insert provided with the kit. Be sure that everything is understood.

The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch and used in the frame provided.

Pipetting of samples and reagents must be done as quickly as possible and in the same sequence for each step.

Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur.

Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells.

Do not let wells dry during assay; add reagents immediately after completing the rinsing steps. Allow the reagents to reach room temperature (21-26°C) before starting the test. Temperature will affect the absorbance readings of the assay. However, values for the patient samples will not be affected.

Never pipet by mouth and avoid contact of reagents and specimens with skin and mucous membranes.

Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may give false results.

Handling should be done in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.

Do not use reagents beyond expiry date as shown on the kit labels.

All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microtiterplate readers.

Do not mix or use components from kits with different lot numbers. It is advised not to exchange wells of different plates even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the plates may result slightly different.

Avoid contact with *Stop Solution* containing 0.5 M H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns.

Some reagents contain Proclin 300, BND and/or MIT as preservatives. In case of contact with eyes or skin, flush immediately with water.

TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.

Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guideline or regulation.

For information on hazardous substances included in the kit please refer to Material Safety Data Sheets. Material Safety Data Sheets for this product are available upon request directly from Biovendor.

## 6. TECHNICAL HINTS

---

All reagents and specimens must be allowed to come to room temperature before use.

All reagents must be mixed without foaming.

Once the test has been started, all steps should be completed without interruption.

Use new disposal plastic pipette tips for each calibrator, control or sample in order to avoid cross contamination

Absorbance is a function of the incubation time and temperature. Before starting the assay, it is recommended that all reagents are ready, caps removed, all needed wells secured in holder, etc. This will ensure equal elapsed time for each pipetting step without interruption.

As a general rule the enzymatic reaction is linearly proportional to time and temperature.

## 7. REAGENT SUPPLIED



12x8 (break apart) strips, 96 wells

Wells coated with anti- human Placental Lactogen antibody (monoclonal)

CAL	0
-----	---

Zero Calibrator (also used as Sample Diluent), 1 vial, 90 mL, ready to use

0 mg/L

contain non-mercury preservative.

CAL	N
-----	---

4 vials, 0.5 ml, ready to use

Concentrations: 1.25 – 5.0 – 10.0 – 20 mg/L

Conversion: 1 mg/L = 1 mIU/L

*The calibrators are calibrated against the NIBSC International Standard for human Placental Lactogen IRP (73/545) contain non-mercury preservative.*

CONTROL	N
---------	---

2 vials, 0.5 ml, ready to use

For control values and ranges please refer to vial label or Quality Control Sheet.

Contain non-mercury preservative.

Ab	HRP
----	-----

Enzyme Conjugate, 1 vial, 11 mL, ready to use

Anti- human Placental Lactogen antibody conjugated to horseradish peroxidase

contain non-mercury preservative.

CHROM	TMB
-------	-----

Substrate Solution, 1 vial, 14 mL, ready to use

Tetramethylbenzidine (TMB)

STOP	STOLN
------	-------

Stop Solution, 1 vial, 14 mL, ready to use

contains 0.5M H<sub>2</sub>SO<sub>4</sub>,

Avoid contact with the stop solution. It may cause skin irritations and burns.

**Note:** Additional *Zero Calibrator* for sample dilution is available on request.

## 8. MATERIAL REQUIRED BUT NOT SUPPLIED

---

A microtiter plate calibrated reader ( $450\pm 10$  nm).  
Calibrated variable precision micropipettes.  
Absorbent paper.  
Aqua dest.  
Tubes for sample dilution (12x75 mm)

## 9. PREPARATION OF REAGENTS

---

Allow all reagents and required number of strips to reach room temperature prior to use.

### »» Disposal of the Kit

The disposal of the kit must be made according to the national regulations. Special information for this product is given in the Material Safety Data Sheets.

### »» Damaged Test Kits

In case of any severe damage of the test kit or components, Biovendor have to be informed written, latest one week after receiving the kit. Severely damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed according to the official regulations.

## 10. PREPARATION OF SAMPLES

---

Serum must be used in this assay.

Do not use haemolytic, icteric or lipaemic specimens.

Please note: Samples containing sodium azide should not be used in the assay.

- **Specimen collection**

### Serum

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time.

- **Specimen storage**

Specimens should be capped and may be stored for up to 5 days at 2°C to 8°C prior to assaying. Specimens held for a longer time should be frozen only once at -20°C prior to assay. Thawed samples should be inverted several times prior to testing.

- **Specimen dilution**

Before starting the assay the sample must be **prediluted 1:100 with Zero Calibrator**:

1. Add for each sample 1mL of Zero Calibrator in one test tube.
2. Add 10 $\mu$ L of each sample to the appropriate test tube. Mix all tubes for 10 seconds on a Vortex mixer (avoid foaming).

The internal controls (Control Low and High) are ready to use and need not to be diluted. External controls (such as Bio-Rad controls) has to be handled like samples.

If in an initial assay, a specimen is found to contain more than the highest calibrator, the specimens can be further diluted with *Zero Calibrator* and reassayed as described in Assay Procedure.

For the calculation of the concentrations this dilution factor has to be taken into account.

## 11. ASSAY PROCEDURE

---

### NOTE:

Manual Pipetting: It is recommended that no more than 32 wells be used for each assay run. Pipetting of all calibrators, samples, and controls should be completed within 3 minutes.

Automated Pipetting: A full plate of 96 wells may be used in each assay run. However, it is recommended that pipetting of all calibrators, samples, and controls be completed within 3 minutes.

Each run must include a calibration curve.

1. Secure the desired number of Microtiter wells in the holder.
2. Dispense **10  $\mu$ L** of each *Standard, control and prediluted samples* and prediluted samples with new disposable tips into appropriate wells.
3. Dispense **100  $\mu$ L** *Enzyme Conjugate* into each well.
4. Thoroughly mix for 10 seconds. It is important to have a complete mixing in this step.
5. Incubate for **30 minutes** at room temperature.
6. Briskly shake out the contents of the wells.  
Rinse the wells 5 times with distilled water (300  $\mu$ L per well). Strike the wells sharply on absorbent paper to remove residual droplets.

### **Important note:**

The sensitivity and precision of this assay is markedly influenced by the correct performance of the washing procedure!

7. Add **100  $\mu$ L** of *Substrate Solution* to each well.
8. Incubate for **10 minutes** at room temperature.
9. Stop the enzymatic reaction by adding **50  $\mu$ L** of *Stop Solution* to each well.  
Read the OD at **450 $\pm$ 10 nm** with a microtiter plate reader **within 10 minutes** after adding the *Stop Solution*.



## 12. CALCULATIONS

---

1. Calculate the average absorbance values for each set of calibrators, controls and patient samples.
2. Using semi-logarithmic graph paper, construct a calibration curve by plotting the mean absorbance obtained from each calibrator against its concentration with absorbance value on the vertical(Y) axis and concentration on the horizontal (X) axis.
3. Using the mean absorbance value for each sample determine the corresponding concentration from the calibration curve.
4. Automated method: The results in the IFU have been calculated automatically using a 4 PL 4 Parameter curve fit. (4 Parameter Rodbard or 4 Parameter Marquardt are the preferred methods). Other data reduction functions may give slightly different results.
5. The concentration of the samples can be read directly from this calibration curve. Samples with concentrations higher than that of the highest calibrator have to be further diluted or reported as > 20mg/L. For the calculation of the concentrations this dilution factor has to be taken into account.

### Example of typical standard curve

The following data is for demonstration only and **cannot** be used in place of data generations at the time of assay.

Calibrator	Optical Units (450 nm)
Zero Calibrator (0 mg/L)	0.03
Calibrator 1 (1.25 mg/L)	0.17
Calibrator 2 (5.0 mg/L)	0.65
Calibrator 3 (10.0 mg/L)	1.17
Calibrator 4 (20.0 mg/L)	1.84

### Expected values

It is strongly recommended that each laboratory should determine its own normal and abnormal values.

In a study, using the Biovendor human Placental Lactogen ELISA, the following values are observed:

<b>Pregnancy week</b>	<b>Value [mg/L]</b>
10. - 12.	0.05 - 1.00
12. - 14.	0.10 - 1.7
14. - 16.	0.3 - 2.8
16. - 18.	0.5 - 3.5
18. - 20.	0.9 - 4.0
20. - 22.	1.1 - 5.0
22. - 24.	1.3 - 5.8
24. - 26.	1.6 - 6.7
26. - 28	2.0 - 7.7
28. - 30.	2.7 - 8.5
30. - 32.	3.2 - 9.5
32. - 34.	3.7 - 10.1
34. - 36.	4.0 - 10.7
36. - 38.	4.3 - 11.2
38. - 40	4.4 - 11.7
40. - 42.	4.3 - 11.6

**Example:**

If the patient is in the first half of pregnancy week 12 please use the range for week 10-12. If the patient is in the second half of week 12 then the range for weeks 12-14 should be used.

The results alone should not be the only reason for any therapeutic consequences. The results should be correlated to other clinical observations and diagnostic tests.

### 13. PERFORMANCE CHARACTERISTICS

- **Assay Dynamic Range**

The range of the assay is between 0.04 – 20 mg/L.

- **Specificity of Antibodies (Cross Reactivity)**

The following substances were tested for cross reactivity of the assay:

Antigen tested	Equivalent to human Placental Lactogen
hCG 2000IU/L	undetectable
AFP 300KIU/L	undetectable
hGH 100 µg/L	undetectable
Prolactin 200 µg/L	undetectable

- **Analytical Sensitivity**

The analytical sensitivity was calculated from the mean plus two standard deviations of twenty (20) replicate analyses of the *Zero Calibrator* and was found to be 0.043 mg/L.

- **Precision**

*Intra Assay Variation* (n=18)

The within assay variability is shown below:

Sample	Mean (mg/L)	CV (%)
1	0.66	6.06
2	2.34	5.55
3	6.24	6.73

*Inter Assay Variation*

The between assay variability is shown below:

Sample	n	Mean (mg/L)	CV (%)
1	39	0.68	8.82
2	24	2.52	7.14
3	24	6.87	5.67

- Recovery**

Samples have been spiked by adding human Placental Lactogen solutions with known concentrations to three different sera containing different amounts of endogenous analyt. The % Recovery has been calculated by multiplication of the ratio of the measurements and the expected values with 100..

Sample	Endogenous human Placental Lactogen (mg/L)	Added Conc. (mg/L)	Measured Conc. (mg/L)	Expected Conc. (mg/L)	Recovery (%)
1	0.00	0.00	0.00		
		0.63	0.55	0.63	88.60
		2.50	2.40	2.50	96.20
		5.00	5.21	5.00	104.20
		10.00	8.58	10.00	85.80
2	1.93	0.00	1.93		
		0.63	1.39	1.59	87.60
		2.50	3.16	3.46	91.40
		5.00	5.21	5.96	87.40
		10.00	9.81	10.96	89.50
3	4.67	0.00	4.67		
		0.63	2.56	2.96	86.40
		2.50	4.46	4.83	92.30
		5.00	6.69	7.33	91.20
		10.00	11.66	12.33	94.50

- Linearity**

Sample	Dilution	Measured Conc. (mg/L)	Expected Conc. (mg/L)	Recovery (%)
1	None	1.93	1.93	
	1:2	0.84	0.96	86.7
	1:4	0.45	0.48	93.7
	1:8	0.26	0.24	108.0
	1:16	0.12	0.12	95.5
2	None	4.67	4.67	
	1:2	2.28	2.33	97.5
	1:4	1.06	1.17	90.7
	1:8	0.64	0.58	109.2
	1:16	0.31	0.29	106.3

## 14. LIMITATIONS OF USE

---

Reliable and reproducible results will be obtained when the assay procedure is performed with a complete understanding of the package insert instruction and with adherence to good laboratory practice. Any improper handling of samples or modification of this test might influence the results.

- **Interfering Substances**

Haemoglobin (up to 4 mg/mL), Bilirubin (up to 0.5 mg/mL) and Triglyceride (up to 30 mg/mL) have no influence on the assay results.

- **Drug Interferences**

Until today no substances (drugs) are known to us, which have an influence to the measurement of human Placental Lactogen in a sample.

- **High-Dose-Hook Effect**

No hook effect was observed in this test up to 700 mg/L of human Placental Lactogen.

## 15. QUALITY CONTROL

---

- Good laboratory practice requires that controls be run with each calibration curve. A statistically significant number of controls should be assayed to establish mean values and acceptable ranges to assure proper performance.
- It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results. Use controls at both normal and pathological levels.
- The controls and the corresponding results of the QC-Laboratory are stated in the Quality Control Sheet added to the kit. The values and ranges states on the Quality Control Sheet always refer to the current kit lot and should be used for direct comparison of the results.
- It is also recommended to make use of national or international Quality Assessment programs in order to ensure the accuracy of the results.
- Employ appropriate statistical methods for analysing control values and trends. If the results of the assay do not fit to the established acceptable ranges of control materials patient results should be considered invalid.
- In this case, please check the following technical areas: Pipetting and timing devices; photometer, expiration dates of reagents, storage and incubation conditions, aspiration and washing methods.
- After checking the above mentioned items without finding any error contact your distributor or Biovendor directly.

## 16. LEGAL ASPECTS

---

### **Reliability of Results**

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.

The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications. In case of any doubt or concern please contact Biovendo.

### **Therapeutical Consequences**

Therapeutical consequences should never be based on laboratory results alone even if all test results are in agreement with the items as stated under point 11.1. Any laboratory result is only a part of the total clinical picture of a patient.

Only in cases where the laboratory results are in acceptable agreement with the overall clinical picture of the patient should therapeutical consequences be derived.


The test result itself should never be the sole determinant for deriving any therapeutical consequences.


### **Liability**

Any modification of the test kit and/or exchange or mixture of any components of different lots from one test kit to another could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement.

Claims submitted due to customer misinterpretation of laboratory results subject to point 11.2. are also invalid. Regardless, in the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

NOTES

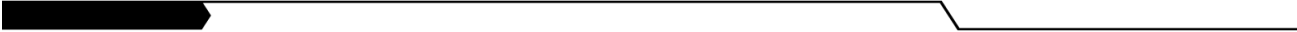

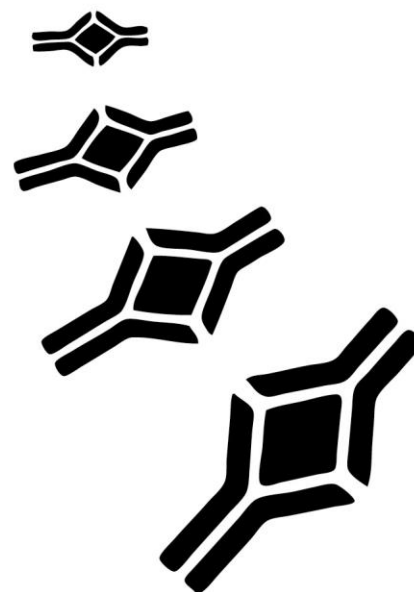

**BioVendor – Laboratorní medicína, a.s.**

Karasek 1767/1, 621 00 Brno, Czech Republic

Phone: +420-549-124-185, Fax: +420-549-211-460

E-mail: [info@biovendor.com](mailto:info@biovendor.com), [sales@biovendor.com](mailto:sales@biovendor.com)

Web: [www.biovendor.com](http://www.biovendor.com)

There are BioVendor branches and distributors near you.

To find the office closest to you, visit [www.biovendor.com/contact](http://www.biovendor.com/contact)

