

## Fibroblast Growth Factor 19 Human E. coli

### Product Data Sheet

**Type:** Recombinant

**Source:** E. coli

**Species:** Human

**Other names:** FGF-19, UNQ334/PRO533

**Cat. No.:**

RD172107100

(0.1 mg)

### Description

Total 206 AA. Mw: 23 kDa (calculated). UniProtKB acc.no. O95750. N-terminal His-tag, 14 extra AA (highlighted).

### Introduction to the Molecule

The FGFs are a family of more than 20 small (~17-26 kDa) secreted peptides. The initial characterization of these proteins focused on their ability to stimulate fibroblast proliferation. This mitogenic activity was mediated through FGF receptors (FGFRs) 1, 2, or 3. A fourth closely related tyrosine kinase receptor (FGFR4) was able to bind the FGFs but did not lead to a mitogenic response. FGFs modulate cellular activity via at least 5 distinct subfamilies of high-affinity FGF receptors (FGFRs): FGFR-1, -2, -3, and -4, all with intrinsic tyrosine kinase activity and, except for FGFR-4, multiple splice isoforms, and FGFR-5, which lacks an intracellular kinase domain. There is growing evidence that FGFRs can be important for regulation of glucose and lipid homeostasis. The overexpression of a dominant negative form of FGFR-1 in beta cells leads to diabetes in mice, which thus implies that proper FGF signaling is required for normal beta cell function and glycemia maintenance. FGFR-2 appears to be a key molecule during pancreatic development. Moreover, FGFR-4 has been implicated in cholesterol metabolism and bile acid synthesis.

FGF-19, has been shown to cause resistance to diet-induced obesity and insulin desensitization and to improve insulin, glucose, and lipid profiles in diabetic rodents. Since these effects, at least in part, are mediated through the observed changes in metabolic rates, FGF-19 can be considered as a regulator of energy expenditure.

FGF-21 is preferentially expressed in liver, but an exact knowledge of FGF-21 bioactivity and its mode of action have been lacking to date. FGF-21 is a potent activator of glucose uptake on adipocytes, protects animals from diet-induced obesity when overexpressed in transgenic mice, and lowers blood glucose and triglyceride levels when therapeutically administered to diabetic rodents.

### Research topic

Energy metabolism and body weight regulation

### Amino Acid Sequence

**MRGSHHHHHH** **GMA**SLAFSDA GPHVHYGWGD PIRLRHLYTS GPHGLSSCFL RIRADGVVDC ARGQSAHSLL EIKAVLRV  
AIKGVHSVRY LCMGADGKMQ GLLQYSEEDC AFEERPDG YNVYRSEKHR LPVSLSSAKQ RQLYKNRGFL PLSHFLPMLP  
MVPEEPEDLR GHLESDFSS PLETDSMDPF GLVTGLEAVR SPSFEK


### Source

E. coli

### Purity

>95%

## SDS-PAGE gel

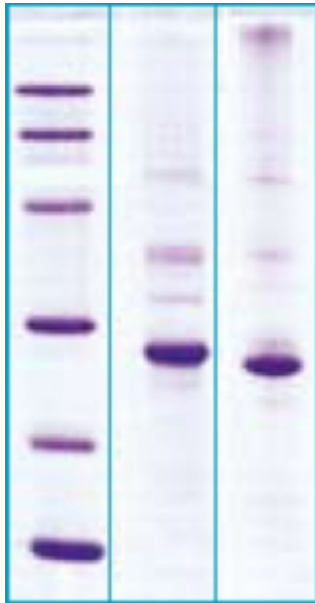

12% SDS-PAGE separation of Human FGF19

1. M.W. marker - 14, 21, 31, 45, 66, 97 kDa

2. reduced and heated sample, 2.5µg/lane

3. non-reduced and non-heated sample, 2.5µg/lane

## Endotoxin

< 1.0 EU/ug

## Formulation

Filtered (0.4 µm) and lyophilized in 0.5 mg/mL in 20mM TRIS, 50mM NaCl, pH 7.5

## Reconstitution

Add deionized water to prepare a working stock solution of approximately 0.5 mg/mL and let the lyophilized pellet dissolve completely. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture.

## Shipping

At ambient temperature. Upon receipt, store the product at the temperature recommended below.

## Storage, Stability/Shelf Life

Store lyophilized protein at -80°C. Lyophilized protein remains stable until the expiry date when stored at -80°C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at -80°C for long term storage. Reconstituted protein can be stored at 4°C for a limited period of time; it does not show any change after one week at 4°C.

## Quality Control Test

BCA to determine quantity of the protein.

SDS PAGE to determine purity of the protein.

LAL to determine quantity of endotoxin.

## Applications

ELISA, Western blotting

## Note

<b>HEADQUARTERS:</b> BioVendor Laboratorní medicína, a.s.	Karasek 1767/1	621 00 Brno CZECH REPUBLIC	Phone: +420-549-124-185 Fax: +420-549-211-460	E-mail: info@biovendor.com sales@biovendor.com Web: www.biovendor.com
AUSTRIA: BioVendor GesmbH	Nußdorfer Straße 20/10	1090 Vienna AUSTRIA	Phone: +43-1-89090-25 Fax: +43-1-89090-2515	E-mail: infoAustria@biovendor.com
GERMANY, SWITZERLAND: BioVendor GmbH	Otto-Hahn-Straße 16	34123 Kassel GERMANY	Phone: +49-6221-433-9100 Fax: +49-6221-433-9111	E-mail: infoEU@biovendor.com
USA, CANADA AND MEXICO: BioVendor LLC	128 Bingham Rd. Suite 1300	Asheville, NC 28806 USA	Phone: +1-828-575-9250 +1-800-404-7807 Fax: +1-828-575-9251	E-mail: infoUSA@biovendor.com
CHINA - Mainland Office: BioVendor Laboratories Ltd	Room 2917, 29/F R & F Ying Feng Plaza No.2 Huaqiang road	Pearl River New Town Guang Zhou, CHINA	Phone: +86-20-38065519 Fax: +86-20-38065529	E-mail: infoCN@biovendor.com