

S100A7A Human E. coli

Product Data Sheet

Type: Recombinant

Source: E. coli

Species: Human

Other names: S100-A15, S100 calcium-binding protein A7A, S100 calcium-binding protein A7-like 1, S100 calcium-binding protein A15

Cat. No.:

RD172213100

(0.1 mg)

Description

Total 109 AA. MW: 12.3 kDa (calculated). UniProtKB acc.no. Q86SG5. N-Terminal His-tag, 9 extra AA (highlighted).

Research topic

Oncology, Others


Amino Acid Sequence

MKHHHHHHAS NTQAERSIIG MIDMFHKYTG RDGKIEKPSL LTMMKENFPN FLSACDKKGI HYLATVFEKK DKNEDKKIDF
SEFLSLLGDI AADYHKQSHG AAPCSGGSQ

Source

E. coli

SDS-PAGE gel

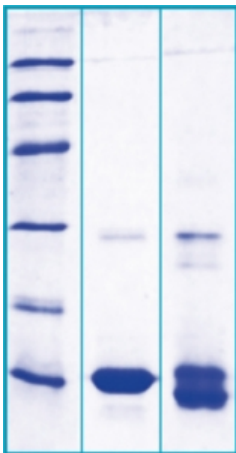

14% SDS-PAGE separation of Human S100A15 protein

1. M.W. marker - 14, 21, 31, 45, 66, 97 kDa

2. reduced and boiled sample, 5µg / lane

3. non-reduced and non-boiled sample, 5µg / lane

Formulation

Filtered (0,4 µm) and lyophilized in 0.5 mg/mL in 20mM TRIS, 50mM NaCl, pH 7.5

Reconstitution

Add deionized water to prepare a working stock solution of approximately 0.5 mg/mL and let the lyophilized pellet dissolve completely. Filter sterilize your culture media/working solutions containing this non-sterile product before using in cell culture.

Shipping

At ambient temperature. Upon receipt, store the product at the temperature recommended below.

Storage, Stability/Shelf Life

Store lyophilized protein at -20°C. Lyophilized protein remains stable until the expiry date when stored at -20°C. Aliquot reconstituted protein to avoid repeated freezing/thawing cycles and store at -80°C for long term storage. Reconstituted protein can be stored at 4°C for a limited period of time; it does not show any change after one week at 4°C.

Quality Control Test

BCA to determine quantity of the protein.

SDS PAGE to determine purity of the protein.

Applications

Western blotting

HEADQUARTERS: BioVendor Laboratorní medicína, a.s.	Karasek 1767/1	621 00 Brno CZECH REPUBLIC	Phone: +420-549-124-185 Fax: +420-549-211-460	E-mail: info@biovendor.com sales@biovendor.com Web: www.biovendor.com
AUSTRIA: BioVendor GesmbH	Nußdorfer Straße 20/10	1090 Vienna AUSTRIA	Phone: +43-1-89090-25 Fax: +43-1-89090-2515	E-mail: infoAustria@biovendor.com
GERMANY, SWITZERLAND: BioVendor GmbH	Otto-Hahn-Straße 16	34123 Kassel GERMANY	Phone: +49-6221-433-9100 Fax: +49-6221-433-9111	E-mail: infoEU@biovendor.com
USA, CANADA AND MEXICO: BioVendor LLC	128 Bingham Rd. Suite 1300	Asheville, NC 28806 USA	Phone: +1-828-575-9250 +1-800-404-7807 Fax: +1-828-575-9251	E-mail: infoUSA@biovendor.com