

BioVendor

Research
and Diagnostic Products



Human Interleukin-36 gamma ELISA

Product Data Sheet

Cat. No.: RAG022R

For Research Use Only

CONTENTS

1.	INTENDED USE	3
2.	HANDLING, STORAGE	3
3.	INTRODUCTION	3
4.	TEST PRINCIPLE	4
5.	TECHNICAL HINTS	4
6.	REAGENT SUPPLIED	5
7.	MATERIALS REQUIRED BUT NOT SUPPLIED	5
8.	PREPARATION OF REAGENTS	5
9.	PREPARATION OF SAMPLES	7
10.	ASSAY PROCEDURE	8
11.	CALCULATIONS	9
12.	PERFORMANCE CHARACTERISTICS	10
13.	TROUBLESHOOTING	11
14.	REFERENCES	12

**»» This kit is manufactured by:
BioVendor – Laboratorní medicína a.s.**

»» Use only the current version of Product Data Sheet enclosed with the kit!

1. INTENDED USE

The Human IL-36 gamma ELISA Kit is to be used for the in vitro quantitative determination of human IL-36 gamma in cell culture supernatants and serum. This ELISA Kit is for research use only.

2. HANDLING, STORAGE

- Reagent must be stored at 2-8°C when not in use.
- Plate and reagents should be at room temperature before use.
- Do not expose reagents to temperatures greater than 25°C.

3. INTRODUCTION

IL-36 α (IL-1F6), IL-36 β (IL-1F8) and IL-36 γ (IL-1F9) are members of the IL-1 cytokine family that bind to IL-36R (IL-1Rrp2) and IL-1RAcP, activating similar intracellular signals as IL-1 and are inhibited by IL-36Ra (1). The expression of IL-36 cytokines occurs mainly in the lung and skin and can be derived from diverse epithelial cell types including keratinocytes, bronchial epithelium as well as macrophages, monocytes and different T cell subsets. IL-36 γ expression is significantly induced upon TLR stimulation and induced in bronchial epithelium exposed to inflammatory stimuli suggesting that these proteins are involved in first-line defences against microorganisms (1,2).

IL-36 family members induce the production of pro-inflammatory cytokines, including IL-12, IL-1 β , IL-6, TNF- α and IL-23, thus promoting neutrophil influx, dendritic cell (DC) activation, polarization of T helper type 1 (Th1) and IL-17-producing T cells ($\alpha\beta$ T cells and $\gamma\delta$ T cells) and keratinocyte proliferation (1). The IL-36 cytokines may represent potential targets for immune-mediated inflammatory conditions or, alternatively, could be used as adjuvants in vaccination. Induction of

IL-36 γ in macrophages upon *M. tuberculosis* infection and its role in the release of antimicrobial peptides has been proposed (2). IL-36 γ is also induced in the lung in various models of asthma and can be produced by bronchial epithelial cells in response to viral infection, smoke or inflammatory cytokines and plays an important role in asthmatic pulmonary inflammation (3).

IL-36 γ might be a potential biomarker of Inflammatory Bowel Disease (4,5), and is highly expressed in inflamed skin from psoriasis patients and in Allergic Contact Dermatitis (6). IL-36 γ serum levels are enhanced and correlate to psoriasis severity and response to treatment with anti-TNF (7).

4. TEST PRINCIPLE

This assay is a sandwich Enzyme Linked-Immunosorbent Assay (ELISA) for quantitative determination of human IL-36 γ in cell culture supernatants and serum. A monoclonal antibody specific for IL-36 γ has been precoated onto the 96-well microtiter plate. Standards (STD) and samples are pipetted into the wells for binding to the coated antibody. After extensive washing to remove unbound compounds, IL-36 γ is recognized by the addition of a biotinylated monoclonal antibody specific for IL-36 γ (DET). After removal of excess biotinylated antibody, streptavidin-peroxidase (STREP-HRP) is added. Following a final washing, peroxidase activity is quantified using the substrate 3,3',5,5'-tetramethylbenzidine (TMB). The intensity of the color reaction is measured at 450 nm after acidification and is directly proportional to the concentration of IL-36 γ in the samples

5. TECHNICAL HINTS

- It is recommended that all standards, controls and samples be run in duplicate.
- Do not combine leftover reagents with those reserved for additional wells.
- Reagents from the kit with a volume less than 100 μ l should be centrifuged.
- Residual wash liquid should be drained from the wells after last wash by tapping the plate on absorbent paper.
- Crystals could appear in the 10X solution due to high salt concentration in the stock solutions. Crystals are readily dissolved at room temperature or at 37°C before dilution of the buffer solutions.
- Once reagents have been added to the 16-well strips, DO NOT let the strips DRY at any time during the assay.
- Keep Substrate Solution protected from light.
- The Stop Solution consists of sulfuric acid. Although diluted, the Stop Solution should be handled with gloves, eye protection and protective clothing.

6. REAGENT SUPPLIED

Kit Components	Quantity
1 plate coated with human IL-36 gamma Antibody	6x16-well strips
2 bottles Wash Buffer 10X	30 ml
2 bottle ELISA Buffer 10X	30 ml
1 vial Detection Antibody	20 μ l
1 vial HRP Labeled Streptavidin (lyophilized)	2 μ g
1 vial human IL-36 gamma Standard (lyophilized)	1 μ g
1 bottle TMB Substrate Solution	12 ml
1 bottle Stop Solution	12 ml
2 plate sealers (plastic film)	
2 silica Gel Minibags	

7. MATERIALS REQUIRED BUT NOT SUPPLIED

- Microtiterplate reader at 450 nm
- Calibrated precision single and multi-channel pipettes. Disposable pipette tips
- Deionized water
- Microtubes or equivalent for preparing dilutions
- Disposable plastic containers for preparing working buffers
- Plate washer: automated or manual
- Glass or plastic tubes for diluting and aliquoting standard

8. PREPARATION OF REAGENTS

NOTE: Prepare just the appropriate amount of the buffers necessary for the assay.

- **Wash Buffer 10X** has to be diluted with deionized water 1:10 before use (e.g. 30 ml Wash Buffer 10X + 270 ml water) to obtain Wash Buffer 1X.
- **ELISA Buffer 10X** has to be diluted with deionized water 1:10 before use (e.g. 10 ml ELISA Buffer 10X + 90 ml water) to obtain ELISA Buffer 1X.

- **Detection Antibody (DET)** has to be diluted to 1:1000 in ELISA Buffer 1X (4 µl DET + 4 ml ELISA Buffer 1X).

NOTE: *The diluted Detection Antibody is not stable and cannot be stored!*

- **HRP Labeled Streptavidin (STREP-HRP)** has to be reconstituted with 100 µl of ELISA Buffer 1X.
 - After reconstitution of STREP-HRP, prepare aliquots and store them at -20°C. Avoid freeze/thaw cycles
 - Dilute the reconstituted STREP-HRP to the working concentration by adding 50 µl in 10 ml of ELISA Buffer 1X (1:200).

NOTE: *The diluted STREP-HRP is not stable and cannot be stored!*

- **Human IL-36 gamma Standard (STD)** has to be reconstituted with 100 µl of ELISA Buffer 1X.
 - This This reconstitution produces a stock solution of 10 µg/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes at room temperature. Mix well prior to making dilutions.

NOTE: *The reconstituted standard is aliquoted and stored at -20°C.*

- Dilute the standard protein concentrate (STD) (10 µg/ml) in ELISA Buffer 1X. A seven-point standard curve using 2-fold serial dilutions in ELISA Buffer 1X is recommended.
- Suggested standard points are:
250, 125, 62.5, 31.2, 15.6, 7.8, 3.9 and 0 pg/ml.

Start with the dilution of the concentrate (STD):

To obtain	Add	Into
100 ng/ml	10µl of IL-36 γ (STD) (10 µg/ml)	990 µl of ELISA Buffer 1X
1 ng/ml	10µl of IL-36 γ (STD) (100 ng/ml)	990 µl of ELISA Buffer 1X

Dilute further for the standard curve:

To obtain	Add	Into
250 pg/ml	250µl of IL-36γ (STD) (1 ng/ml)	750 µl of ELISA Buffer 1X
125 pg/ml	300 µl of IL-36γ (250 pg/ml)	300 µl of ELISA Buffer 1X
62.5 pg/ml	300 µl of IL-36γ (125 pg/ml)	300 µl of ELISA Buffer 1X
31.25 pg/ml	300 µl of IL-36γ (62.5 pg/ml)	300 µl of ELISA Buffer 1X
15.6 pg/ml	300 µl of IL-36γ (31.25 pg/ml)	300 µl of ELISA Buffer 1X
7.8 pg/ml	300 µl of IL-36γ (15.6 pg/ml)	300 µl of ELISA Buffer 1X
3.9 pg/ml	300 µl of IL-36γ (7.8 pg/ml)	300 µl of ELISA Buffer 1X
0 ng/ml	300 µl of ELISA Buffer 1X	Empty tube

9. PREPARATION OF SAMPLES

Serum: Use a serum separator tube. Let samples clot at room temperature for 30 minutes before centrifugation for 20 minutes at 1,000xg. Assay freshly prepared serum or store serum in aliquot at $\leq -20^{\circ}\text{C}$ for later use. Avoid repeated freeze/thaw cycles.

Serum and Cell Culture Supernatant have to be diluted in ELISA Buffer 1X. Samples containing visible precipitates must be clarified before use

NOTE: As a starting point, 1/2 dilution of serum is recommended! If sample values fall outside the detection range of the assay, a lower or higher dilution may be required!


10. ASSAY PROCEDURE

<input type="checkbox"/>	<p>1. Determine the number of 16-well strips needed for the assay and insert them in the frame for current use. The extra strips are left in the bag with 2 silica gel minibags and stored at 4°C.</p> <p>NOTE: <i>Remaining 16-well strips coated with IL-36γ antibody when opened can be stored in the presence of 2 silica gel minibags at 4°C for up to 1 month.</i></p>
<input type="checkbox"/>	<p>2. Add 100 μl of the different standards into the appropriate wells in duplicate! At the same time, add 100 μl of diluted serum or cell culture supernatant samples in duplicate to the wells (see 8.1. Preparation and Storage of Reagents and 8.2 Preparation of Samples).</p>
<input type="checkbox"/>	<p>3. Cover the plate with plastic film and incubate for 2 h at RT°C.</p>
<input type="checkbox"/>	<p>4. Aspirate the coated wells and add 300 μl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of five washes. After the last wash, complete removal of liquid is essential for good performance.</p>
<input type="checkbox"/>	<p>5. Add 100 μl to each well of the diluted Detection Antibody (DET) (see 8.1 Preparation and Storage of Reagents).</p>
<input type="checkbox"/>	<p>6. Cover the plate with plastic film and incubate for 1 h at RT°C.</p>
<input type="checkbox"/>	<p>7. Aspirate the coated wells and add 300 μl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of five washes. After the last wash, complete removal of liquid is essential for good performance.</p>
<input type="checkbox"/>	<p>8. Add 100 μl to each well of the diluted HRP Labeled Streptavidin (STREP-HRP) (see 8.1. Preparation and Storage of Reagents).</p>
<input type="checkbox"/>	<p>9. Cover the plate with plastic film and incubate for 30 min at RT.</p>
<input type="checkbox"/>	<p>10. Aspirate the coated wells and add 300 μl of Wash Buffer 1X using a multichannel pipette or auto-washer. Repeat the process for a total of five washes. After the last wash, complete removal of liquid is essential for good performance.</p>
<input type="checkbox"/>	<p>11. Add 100 μl to each well of TMB substrate solution (TMB).</p>
<input type="checkbox"/>	<p>12. Allow the color reaction to develop at RT in the dark for 15-30 minutes. Do not cover the plate.</p>
<input type="checkbox"/>	<p>13. Stop the reaction by adding 50 μl of Stop Solution (STOP). Tap the plate gently to ensure thorough mixing. The substrate reaction yields a blue solution that turns yellow when Stop Solution (STOP) is added.</p>
	<p>! CAUTION: CORROSIVE SOLUTION !</p>
<input type="checkbox"/>	<p>14. Measure the OD at 450 nm in an ELISA reader.</p>

11. CALCULATIONS

- Average the duplicate readings for each standard, control and sample and subtract the average blank value (obtained with the 0 ng/ml point).
- Generate the standard curve by plotting the average absorbance obtained for each standard concentration on the horizontal (X) axis vs. the corresponding IL-36 γ concentration (pg/ml) on the vertical axis (see 10. TYPICAL DATA).
- Calculate the IL-36 γ concentrations of samples by interpolation of the regression curve formula as shown above in a form of a quadratic equation
- If the test sample was diluted, multiply the interpolated value by the dilution factor to calculate the concentration of human IL-36 γ in the sample.

The following data are obtained using the different concentrations of standard as described in this protocol:

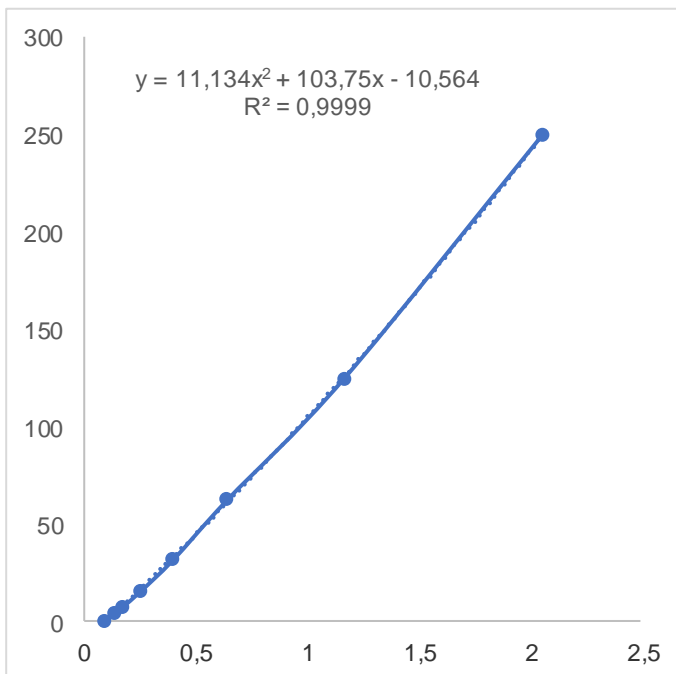

Standard IL-36 γ (pg/ml)	Optical Density (mean)
250	2.0565
125	1.1675
62.5	0.642
31.75	0.3995
15.9	0.2565
7.8	0.176
3.9	0.135
0	0.0955

Figure: Standard curve

12. PERFORMANCE CHARACTERISTICS

- **Sensitivity** (Limit of detection)

The lowest level of IL-36 γ that can be detected by this assay is 3 pg/ml. **NOTE:** The Limit of detection was measured by adding three standard deviations to the mean value of 50 zero standard.

- **Assay Range**

3.9 pg/ml – 250 pg/ml

- **Specificity**

This ELISA is specific for the measurement of natural and recombinant human IL-36 γ . It does not cross-react with human IL-36 α and IL-36 β .

- **Precision:**

Intra-assay:

Four samples of known concentrations of human IL-36 γ were assayed in replicates 8 times to test precision within an assay.

Samples	Means (pg/ml)	SD	CV (%)	n
A1	135.42	6.83	5.05	8
A2	64.28	2.50	3.89	8
A3	16.12	0.72	4.48	8
A4	285.22	11.75	4.12	8

Inter-assay

Three samples of known concentrations of human IL-36 γ were assayed in 4 separate assays to test precision between assays.

Samples	Means (pg/ml)	SD	CV (%)	n
B1	127.23	1.87	1.47	4
B2	29.86	1.01	3.40	4
B3	102.34	8.80	8.60	4
B4	83.78	7.74	9.25	4

- **Linearity**

Different samples containing human IL-36 γ were diluted several fold (1/2 to 1/8) and the measured recoveries ranged from 95% to 106%.

- **Expected values:**

Human IL-36 γ levels range in serum from < 3pg/ml (not detectable) (healthy patient) to > 100pg/ml (psoriasis patient).

13. TROUBLESHOOTING


PROBLEM	POSSIBLE CAUSES	SOLUTIONS
No signal or weak signal	Omission of key reagent	Check that all reagents have been added in the correct order.
	Washes too stringent	Use an automated plate washer if possible.
	Incubation times inadequate	Incubation times should be followed as indicated in the manual.
	Plate reader settings not optimal	Verify the wavelength and filter setting in the plate reader.
	Incorrect assay temperature	Use recommended incubation temperature. Bring substrates to room temperature before use.
High background	Concentration of STREP-HRP too high	Use recommended dilution factor.
	Inadequate washing	Ensure all wells are filling wash buffer and are aspirated completely.
Poor standard curve	Wells not completely aspirated	Completely aspirate wells between steps.
	Reagents poorly mixed	Be sure that reagents are thoroughly mixed.
Unexpected results	Omission of reagents	Be sure that reagents were prepared correctly and added in the correct order.
	Dilution error	Check pipetting technique and double-check calculations.

14. REFERENCES

1. The emergence of the IL-36 cytokine family as novel targets for inflammatory diseases: P.T. Walsh & P.G. Fallon; *Ann. N. Y. Acad. Sci.* (Epub ahead of print) (2016)
2. Role of Interleukin 36 γ in Host Defense Against Tuberculosis: F. Ahsan, et al.; *J. Infect. Dis.* 214, 464 (2016)
3. Rhinovirus-induced modulation of gene expression in bronchial epithelial cells from subjects with asthma: Y.A. Bochkov, et al.; *Mucosal Immunol.* 3, 69 (2010)
4. IL-36 Receptor Promotes Resolution of Intestinal Damage: O. Medina-Contreras, et al.; *J. Immunol.* 196, 34 (2016)
5. Increased Expression of Interleukin-36, a Member of the Interleukin-1 Cytokine Family, in Inflammatory Bowel Disease: A. Nishida, et al.; *Inflamm. Bowel Dis.* 22, 303 (2016)
6. IL-36 γ Is Involved in Psoriasis and Allergic Contact Dermatitis: A. Balato, et al.; *J. Invest. Dermatol.* 136, 1520 (2016)
7. IL-36 γ (IL-1F9) is a biomarker for psoriasis skin lesions: A.M. D'Erme, et al.; *J. Invest. Dermatol.* 135, 1025 (2015)

NOTES



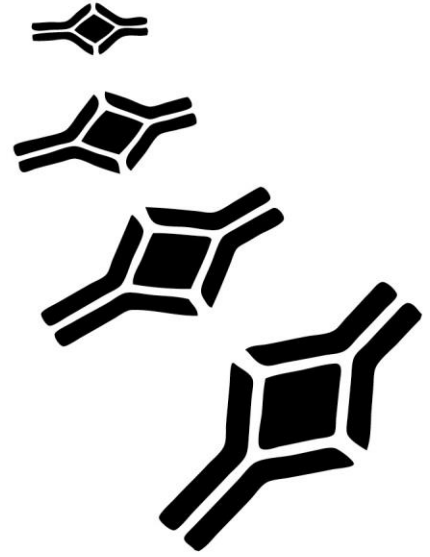

BioVendor – Laboratorní medicína a.s.

Karasek 1767/1, 621 00 Brno, Czech Republic

Phone: +420-549-124-185, Fax: +420-549-211-460

E-mail: info@biovendor.com, sales@biovendor.com

Web: www.biovendor.com

There are BioVendor branches and distributors near you.

To find the office closest to you, visit www.biovendor.com/contact

